




➤ Life sciences & Chemical Analysis



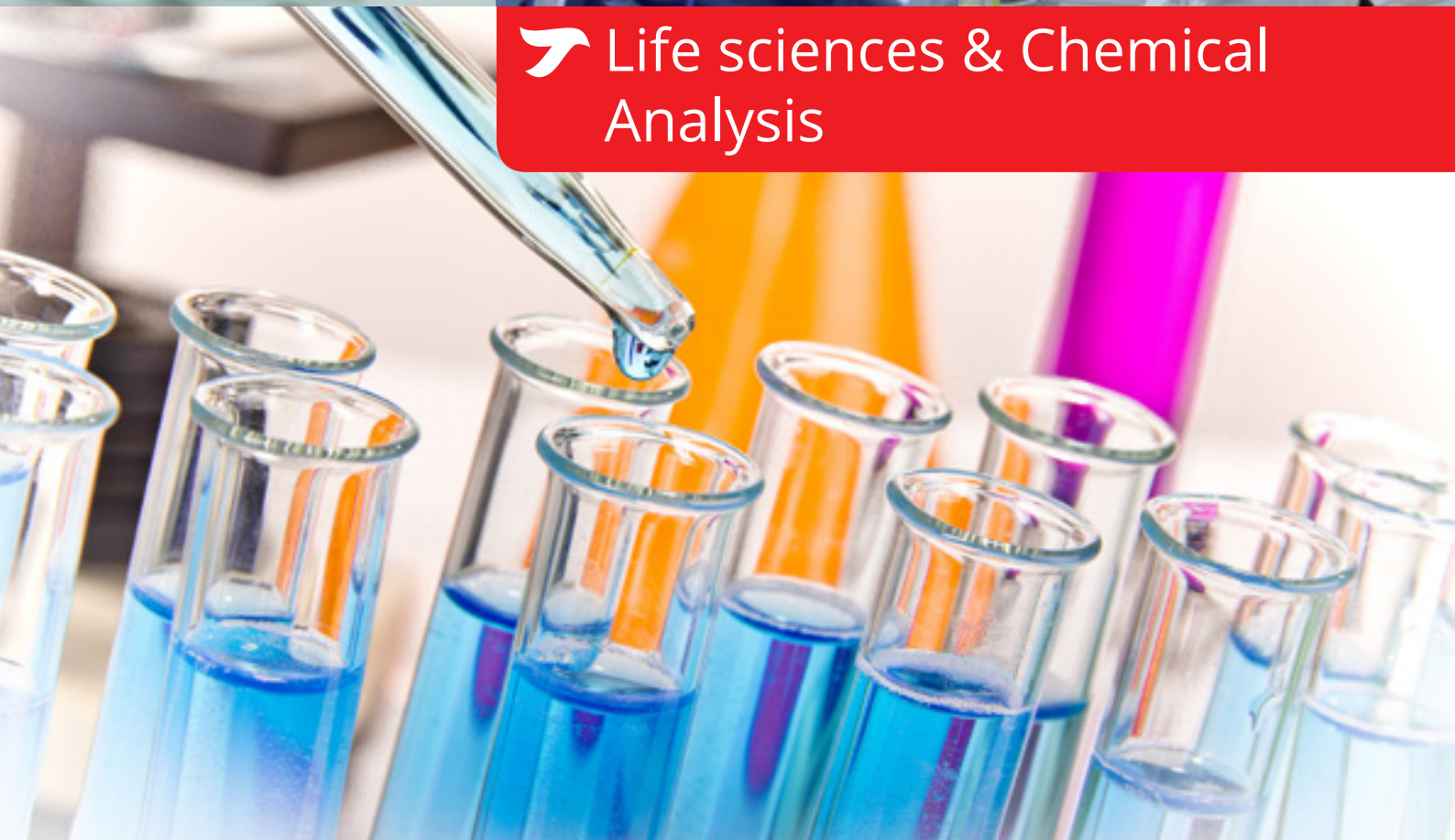

THERMOCOAX

ANALYSIS

IN-house production of MI cable



Applications

- ▶ Temperature Vaporizer inlet
- ▶ Direct inlet probe systems
- ▶ Flexible GC/MS transfer line
- ▶ Adsorption trap
- ▶ Desolvation nozzle
- ▶ Atmospheric pressure chemical
- ▶ Ionization (APCI) probe
- ▶ Ion sources heated probe
- ▶ Differential Scanning Calorimeter (DSC)
- ▶ Differential Thermal Analysis (DTA)
- ▶ Mass Spectrometer probe tip assembly
- ▶ Moisture Sorption analyzers heater
- ▶ Thermal Gravimetric analyzer
- ▶ Gas chromatograph
- ▶ Injectors LC/MS transfer lines
- ▶ Laboratory chemical reactors
- ▶ Reaction vessels sensors

Our capabilities

Our mission is to develop state-of-the-art, complete thermal solutions and this through advanced innovation from conception, to prototype development right the way through to production. We offer tailor-made thermal systems, providing our customers with safer, cleaner and the most cost effective solutions :

- Design / thermal modeling
 - Prototypes
 - Production manufacturing
 - Product support & engineering
 - From cryogenic temperatures up to 1100°C
 - Suitable with atmosphere to UHV environment
 - Thermal modeling / simulation
 - Experience in common and special materials (stainless steels, Inconel® alloys, nickel, tungsten, Aluminum, titanium...)
 - From small up to large size devices
 - Research & testing services
-
- True cold end heater (no connector on the cold/hot junction)
 - Minimum bending radius 3 times the outer diameter
 - High thermal flux > 6 W/cm²
 - High dielectric withstand with an insulation resistance 1.13Ω.m at room temperature
 - Small size heaters 0.5 mm outer diameter available

The industry's
most trusted
Mineral
Insulated
Cable



Separation Technology

Gas Chromatograph (GC)
Liquid Chromatograph (LC)
Ion Chromatograph
Electrophoresis
instrumentation

- ▶ Lower thermal masses with heater OD as low as 0.5 mm
- ▶ Cleanness analytical standards with vacuum brazing technology
- ▶ Custom packaging available with plastic tray and vacuum sealed bags

Gas chromatography PTV inlet



The essentials required for the method are an injection port through which samples are loaded, a «column» on which the components are separated. The temperature of the injection port, column, and detector are controlled by thermostated heaters.

▶ Temperature range


From room temperature to operation at 350°C

▶ Heater

True cold end heaters with integrated thermocouple

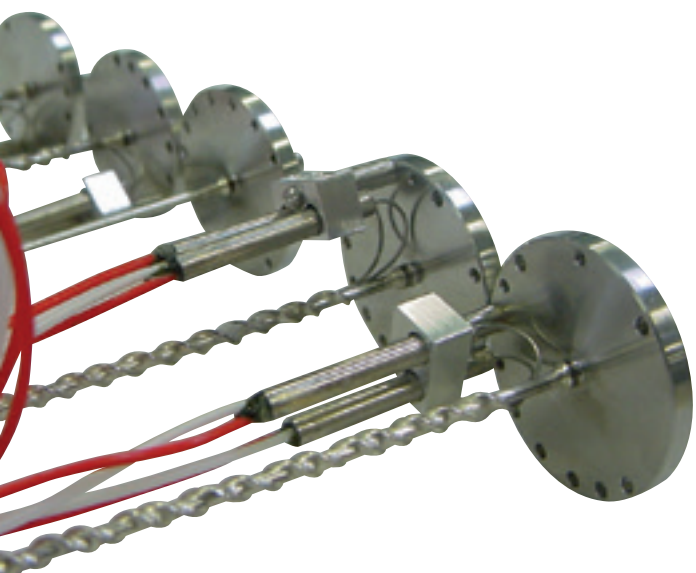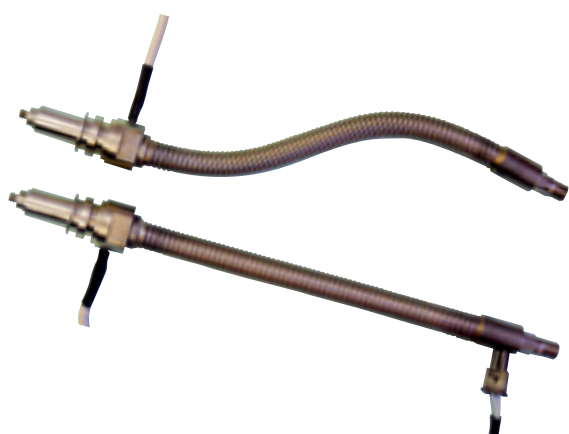

EASE OF OPERATION

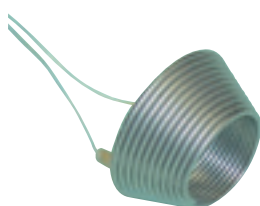

Molecular Analysis


Mass Spectrometer (MS)

- ▣ Quadrupole
- ▣ Ion trap
- ▣ Time of flight
- ▣ Magnetic

Nuclear Magnetic Resonance (NMR)

Moisture sorption analysis

- ▣ Improved heat distribution with customized thermal profile to compensate the heat effect due to the fitting on the machine
- ▣ Design engineered for low voltage power supply
- ▣ 2 wires heater hand wrapped with a different pitch to ensure temperature uniformity

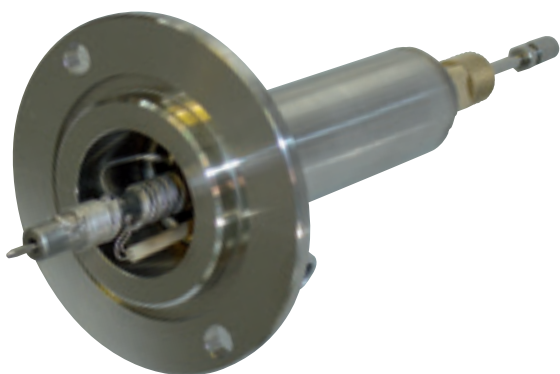

Direct inlet probe

Transfer of low volatile substances into the high vacuum ion source of the MS by a push rod

- ▣ Temperature range

Room temperature to 650°C – 12 V / 5 W Twin core heater with integrated thermocouple type K.

Mass Spectrometer Transfer Line

Rigid and flexible versions for transferring the analytes from the separation unit (GC or LC) into the analytical spectrometer.

MATERIAL

Stainless steel, vacuum brazed heater with varying distribution to obtain temperature uniformity $\pm 2^{\circ}\text{C}$ along full length.

OPERATIONAL VOLTAGE AND POWER


24V – 110W

TEMPERATURE RANGE


Standard operation 200°C - 800°C

HEATER

Special stainless steel or Inconel alloy. The heater is brazed to improve the thermal transfer between the heater and the tube.

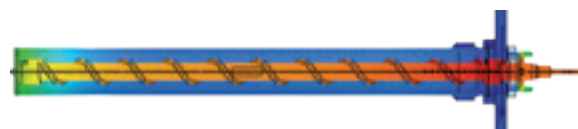

INCREASE RATIO
PERFORMANCE / COST

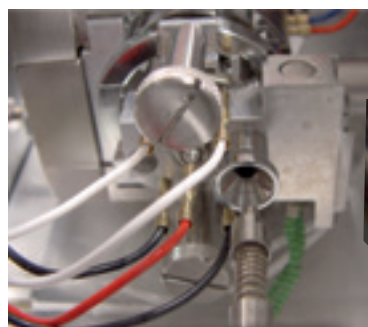

Courtesy of SIM Scientific
Instruments Manufacturer





Thermal Analysis

Differential Scanning Calorimetry (DSC)
Differential Thermal Analysis (DTA)
Determination of thermal conductivity
Thermo gravimetric analysis



- ▶ In house testing with an accredited COFRAC laboratory
- ▶ Complete thermal system design including thermal modelisation, heater wrapped with custom shape and thermocouples for temperature regulation
- ▶ Custom design for faster ramp rates and excellent repeatability



Differential Scanning calorimetry

The sample is heated in an oven at a given power rate until a phase transition is reached.

The power factor used and the time achieved are recorded making the basis of the analytical process. Additional sensors might be added to record other properties during the process.

▶ **Temperature range**

From room temperature up to 850°C

▶ **Heater**

Special stainless steel or Inconel® alloy.

The heater is brazed to obtain the high level of quality and repeatability required by the analytical industry.

SECURE DATE DURING ANALYSIS



Elemental Analysis

CHNOS analysis

Environmental analysis

Agrarian Analysis

MATERIAL

Copper or stainless steel tube with heater coiled and torch brazed along the full length.

OPERATIONAL VOLTAGE/POWER

30 V / 130 W

42 V / 130 W

HEATER ELEMENT


Standard stainless steel sheathed single core heater, Balco® type without cold ends Type K THERMOCOAX thermocouple for regulation



CHNOS Analysis : Adsorption trap

Enrichment of specific gaseous components by trapping on an adsorbate ; followed by thermal desorption for CHNOS analysis.

THERMOCOAX heaters provide temperature uniformity in the dynamic separation system while the gas are retarded in the specific adsorption traps.

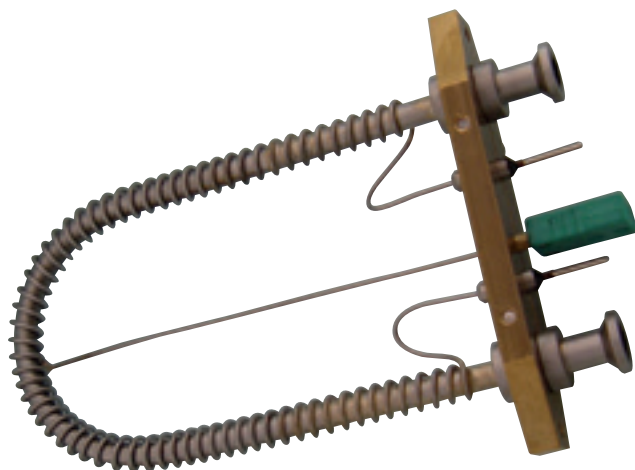

Special adsorption traps extend the dynamic range of the vario EL III (Elementar Analysensysteme GmbH)



EXPAND END USER APPLICATION BASE

Concept, design and computerizing model techniques (CMT) to simulate the prototype working performances

Inhouse testing and calibration using COFRAC certified laboratory

R&D department for static and dynamic thermal modeling



OUR COMPANY

With nearly 60 years of experience in heating solutions and temperature measurement, THERMOCOAX has acquired a great deal of skill and expertise.

THERMOCOAX products are widely used and endorsed in many industries where the highest quality and utmost reliability are essential. All our mineral insulated cables are manufactured in-house with our proprietary and unique procedures.



sales@thermocoax-analysis.com
www.thermocoax-analysis.com

YOUR CONTACT

FRANCE :

THERMOCOAX SAS
 40 Bd Henri Sellier
 F 92156 SURESNES Cedex
 Tel.: +33 1 41 38 80 50
 Fax: +33 1 41 38 80 70
info@thermocoax.com

GERMANY

THERMOCOAX ISOPAD GmbH
 Englerstrasse 11
 D-69 126 HEIDELBERG
 Tel.: +49 6221 3043-0
 Fax: +49 6221 3043-956
Isopad.info@thermocoax.com

USA

THERMOCOAX Inc.
 6825 Shiloh Road East,
 Ste B-3
 ALPHARETTA, GA 30005
 Tel.: +1 678 947 5510
 Fax: +1 678 947 4450
info@thermocoax.us

UK

Tel.: +44 (0)121 767 1822
info-uk@thermocoax.com

CHINA

Tel.: 13701325459
info-china@thermocoax.com



Site Planquignon



Thermocouples
 & Heating Elements
 ATEX
 certified and notified

